

ClearView™

Uterine manipulator injector

The ClearView™ provides full ante version, retroversion and lateral motion for greater visibility and movement.

The device contains a pivot at the cervical OS for focused control. The ClearView™ allows for easy access and simultaneous use with laparoscope without breaking the sterile barrier. It can also provide full exposure for better vision which allows for quicker and safer procedures.

ClearView™ is both a uterine manipulator and a uterine injector designed for single use. The Clearview™ is made with an elastomeric tip to reduce perforation risk and can perform all manipulations with patient in the supine position. No tenaculum is required in most cases and the device results in a reduced pressure on adnexal tissue. A narrow tip reduces the need for dilation and decreases tissue trauma. With the ClearView™ no assistant is needed to hold the manipulator. The device has a balloon seal for chromotubation.

The extended length of the ClearView™ is 345mm (13.6 inches) and the outer diameter of the distal tip is 4.5mm (0.18 inches).

Technical Specifications

This product is designed as a single unit body with a rotation control knob which controls the tip. It has:


- an enclosure, a dull lumen obturating tip which pivots at the cervical OS;
- a soft inflatable uterine cavity balloon with external inflation tubing for securing the device in the uterine cavity;
- a luer fitting to accommodate a syringe for contrast media injection through the distal tip assembly;

Ordering Information

| | |
|-------------|--------------------------------|
| Part Number | CV-1000 (7cm tip) |
| | CV-2000 (9cm tip) |
| Quantity | 5 units per case |
| Description | Clearview™ Uterine Manipulator |
| Delivery | Usually ex-stock |

ColpoCups (CVC-2000)

| |
|-------------------------------------|
| For use with CV1000, Cv2000 in LAVH |
| Please call for availability |

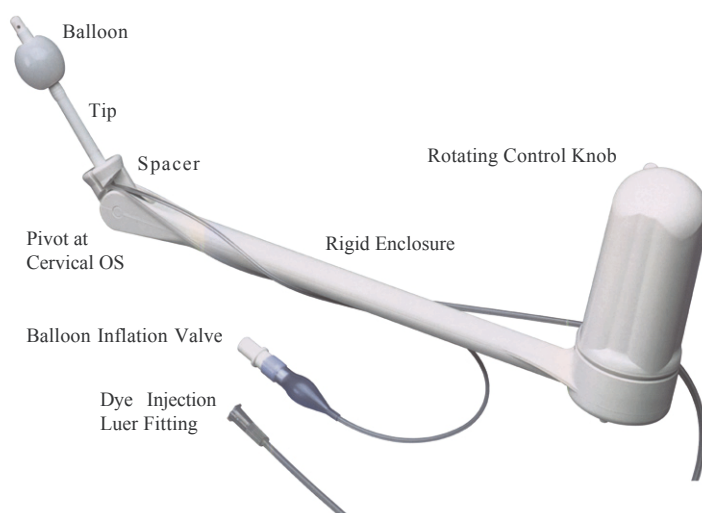

www.clinicalinnovationseurope.com

Distributed