


Catheters for sonohysterography, hysterosalpingography

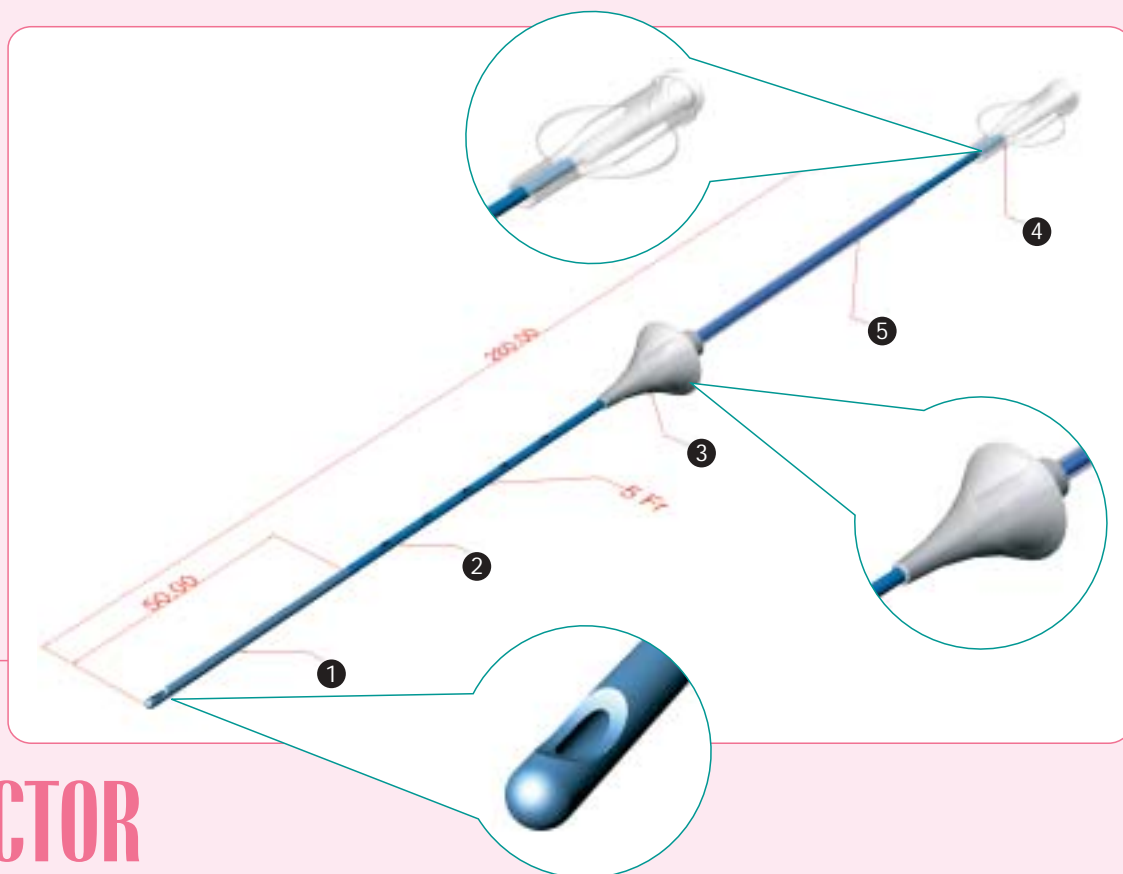

MINIJECTOR

Iniettore uterino per sonoisterografia
Sonohysterography uterine catheter

● Per l'iniezione di soluzione salina nelle procedure di sonoisterografia.

Il dispositivo è costituito da:

- un catetere azzurro centimetrato sulla parte distale (2), per un posizionamento facilitato;
- una punta radiopaca morbida con foro all'estremità (1) per permettere l'introduzione della soluzione salina;
- un introduttore rigido blu (5) che facilita l'inserimento nel collo dell'utero dello stopper costruito in morbido materiale plastico (3) che ne assicura la perfetta tenuta.

Il catetere termina con un raccordo luer lock femmina (4) per la connessione della siringa.

● To infuse saline solution into the uterine cavity in sonohysterography diagnostic procedure. The device is composed by:

- a light blue catheter with centimeter markings on the distal tip (2) for easy positioning into the uterine cavity.
- a radiopaque side-holed soft tip (1) to allow saline infusion.
- a dark blue rigid introducer (5) to advance soft rubber stopper (3) into uterus neck for a guarantee of perfect sealing. The catheter ends in a female luer-lock connector (4) for syringe connection.

Cod.	Descrizione - Description
770314	Ø esterno 1,67 mm. - O.D. 1,67 mm. (5 Fr)

Scatola da 15 pezzi. Fornito in confezione singola, monouso, sterilizzato a EtO e marcato CE.

15 pcs. per box. Single packaged, single use, EtO sterilized and CE marked.